

# Earol Swim<sup>®</sup> Tea Tree Oil

Earol Swim<sup>®</sup> Tea Tree Oil delivers a metered dose of a unique blend of Olive Oil & Tea Tree Oil into the outer ear finely coating the auditory canal creating a water resistant barrier.

Olive Oil & Tea Tree Oil have natural anti-fungal, anti-bacterial, anti-viral, and anti-infection properties and when combined with the pump action metered spray delivery system will help prevent the incidence of Trapped Water & Swimmer's Ear.

Please read the directions for use carefully before you use Earol Swim<sup>®</sup> Tea Tree Oil.

## **BEFORE USING A NEW BOTTLE OF EAROL SWIM<sup>®</sup> TEA TREE OIL.**

■ Before using a new bottle of Earol Swim Tea Tree Oil, warm in palm of hand, shake well, then press the actuator down 3 or 4 times until you get a spray (See Diagram b). When doing this do not point the actuator at yourself or anyone else. This procedure will prime the system and ensure the product is ready to use.

## **DIRECTIONS FOR USE**

- Apply Earol Swim<sup>®</sup> Tea Tree Oil to both ears before entering the water.
- Pull the upper ear or Pinna upwards and backwards (see diagram a)
- Holding the bottle in an upright position with the nozzle placed into the ear canal press the actuator firmly. (A single actuation per ear should be sufficient.)
- The actuator button (see diagram b) can be removed and rinsed in hot water after use and then seated back on the stem.

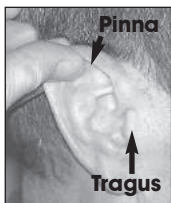


Diagram a



Diagram b

# Earol Swim<sup>®</sup> Tea Tree Oil

## **ADDITIONAL INFORMATION**

- Olive Oil is a natural product and has wide daily applications, as such no side-effects have been reported with this product or with this product interacting with other medicines. However this product should not be used if you have a sensitivity or allergy to Olive Oil.
- Tea Tree oil is a natural product and has been shown to have anti-septic and anti-bacterial properties. No side-effects have been reported with this product when used in the percentages formulated within this product or with this product interacting with other medicines. However this product should not be used if you have a sensitivity or allergy to Tea Tree Oil.
- There are no known harmful side-effects when using this product during pregnancy and/or breast feeding, however you should always inform your Doctor before use.
- There are no known side effects from over or under dosing with this product.
- If you feel you have experienced an allergic reaction stop using this product and inform your Doctor or pharmacy immediately.

## **WARNINGS**

- Earol Swim<sup>®</sup> Tea Tree Oil is for single person use only.
- If you have or have had pre-existing outer ear conditions, it is recommended that you speak with you Doctor or practice Nurse prior to using this product.
- Always store Earol Swim<sup>®</sup> Tea Tree Oil spray in the case provided at normal room temperature.
- The product should be discarded six months after first application.
- Dispose of the product as per Local Authority advice.
- Do not spray near naked flame.
- Earol Swim<sup>®</sup> Tea Tree Oil is non-medicinal, however always remember to keep all medicines in a safe place where children cannot reach them.
- Earol Swim<sup>®</sup> Tea Tree Oil is only for use in the ear canal. Do not spray near eyes, nose or mouth.

**Earol Swim<sup>®</sup> Tea Tree Oil is manufactured in the UK by  
HL Healthcare Ltd, 10 Hazelwood Road, Hale, Cheshire WA15 9AX.**  
Earol Swim<sup>®</sup> is a registered trademark of HL Healthcare.  
UK and International Patent: 1778315